




Digital technology improves life

# P53

Smarter, easier, more effective ways  
to partner with your lab

## Smart Dental Milling Machines


Automatic Calibration And Great Durability

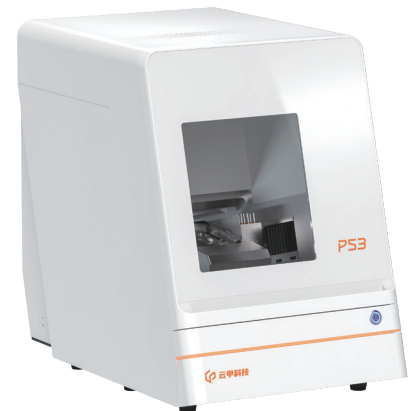



## Automatic Calibration

Calibration with one click. Automatic calibration ensures the milling machine in high-precision milling mode.

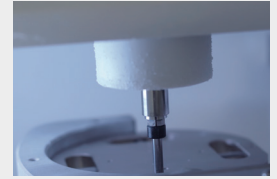

## Save Materials

New C-clamp design to improve the material utilization rate at the edge.

Crown, bridge, coping, implant, full denture and other restoration cases. Support multiple materials, zirconia, wax, PEEK, PMMA, etc.



## 01 Supports the 90° Vertical Milling

The 90° vertical milling method can make the design of the front teeth more aesthetic.

## 02 Independent R&D

All hardware and software are developed by Up3d. Intelligent control system perfectly merged with durable and integrated design.

## 03 Five Smart Points

A smart dental milling machine that defines the era of intelligent dental laboratories.


## 04 UPCNC Software

One interface for multi-mills, easy to operate and one-click milling.

# Maximizing your efficiency

## P53- 5-axis milling machine

A 5-axis milling machine developed by UP3D, equipped with intelligent features such as automatic calibration, monitoring of burs life, and automatic renewal, the P53 digitizes and scales up traditional processing methods.

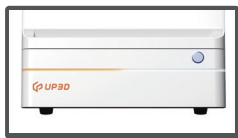

### Led Light Milling Progress Bar

Displaying the milling progress of P53 can reduce personnel inspection time and improve efficiency.



### Remote Connection

One software controls all milling machines, and enables real-time monitoring of the production status of multiple machines.



### Tools Lifetime Display

P53 is able to realize automatic calculation of tools lifetime, reasonable planning of tools life cycle, and reduce the risk of zirconia breakage and milling failure with UPCAM.



### Resume

Automatically continue milling from the last interrupted step, avoiding repeated work and improving work efficiency by 80%.

# 90° Vertical Milling

## WHAT IS THE 90° MILLING & DESIGN?

Ensuring superior results when milling interproximal areas, anterior teeth, and gingival lines




### YOUR SATISFACTION IS OUR PURSUIT

Quality assurance lab testing all products before delivery.




### PROFESSIONAL TECHNICAL SUPPORT

24-hour online technical support and service to meet your needs

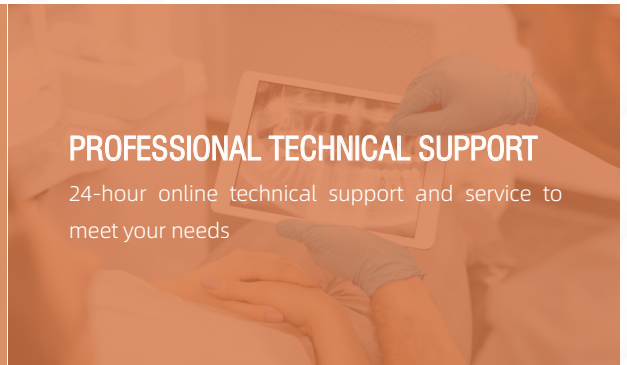


### WARRANTY PERIOD

One year for hardware, everlasting for software

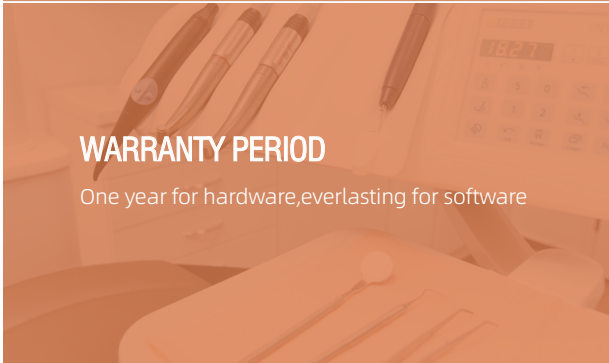

### AFTER-SALES GUARANTEE


No annual fee for software, and professional technical team do the excellent service.

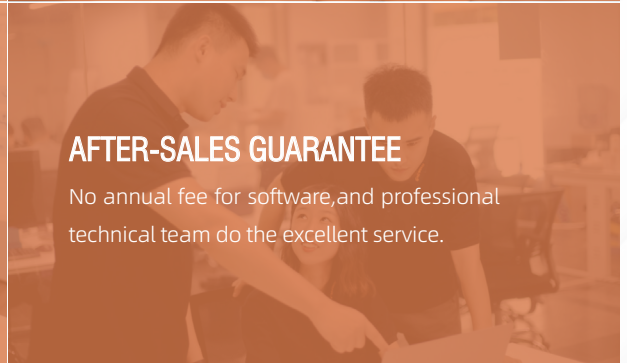

Digital technology improves life



## UP3D 5-Axis Mill Machine

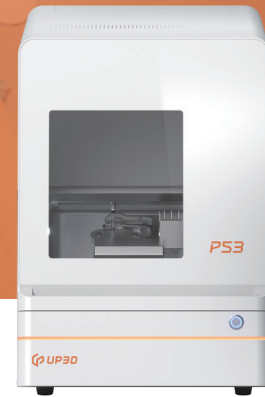

Model	P53
Dimension	443.5x718x628.5 mm
Weight	87kg
Input voltage	AC 200~240V 50~60Hz
Maximum power	1.1KW
Spindle power	0.35KW at max
Axis quantity	5-axis simultaneous motion
Milling range	X/Y/Z:167/206/104mm A:±30° B:±360°
Milling mode	Dry milling
Max rotation speed	30,000rpm
Max feed rate	3000mm/min
Tools quantity	11
Tools changing method	Automatic (Air pressure >0.5MPa )
Spindle cooling	Air cooling with compression air (Air pressure 0.25~0.35MPa)
Milling materials	Zirconia, Wax, PMMA, PEEK and soft materials
Milling time	Coping : Zirconia-10min Wax- 4.5min Crown : Zirconia-13min Wax-5.5min
Tools diameter	2 + 1 + 0.6 ( Φ4 calibration bur )
Disc specification	Round disc with Φ98, 30mm max for thickness
Spindle clamping diameter	Φ4.0mm
CAM supported	UPCAM and Millbox
Data transmission	Network interface





UP3D is a high tech company and the few companies in the world that are capable of developing a full dental CAD/CAM solution including software and hardware, from scanning, designing, to milling.

UP3D remains dedicated to our customers' requirements for training and support as well as for innovative, high-quality products that perform reliably. We have native technical support to deal with any after-sales problems for you to ensure that you can use our products happily.

## CONTACT US

 [www.up3ds.com](http://www.up3ds.com)

 [info@up3d.cn](mailto:info@up3d.cn)

 +86-755-26983202

